



Aegis® Product Line with ACP


Aegis products contain the “smart material” Amorphous Calcium Phosphate (ACP). ACP slowly releases calcium and phosphate ions, the basic building blocks of teeth. ACP is referred to as a “smart material” because it only releases calcium and phosphate ions when the surrounding pH drops to a level where it could start to dissolve the tooth. Once released, the ACP will act to neutralize the acid and buffer the pH. ACP acts as a reinforcement to the tooth’s natural defense system only when it’s needed.



Aegis® Crown & Bridge Cement with ACP



A syringe dispensed, two-paste, dual-cure modified resin polymer material. It is the first crown and bridge resin cement containing an Amorphous Calcium Phosphate (ACP) filler. ACP has demonstrated the ability to reduce micro-leakage related decay due to its ability to accelerate the tooth’s own remineralization process. Ideal for use as a permanent luting cement for crown and bridge, inlay, and veneer cementation applications. If exposed to a curing light, Aegis C&B cures immediately and has post-setting chemical curing to provide improved properties. New activator makes this the perfect product!

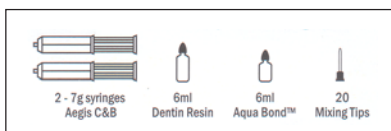

Standard Kit
Shade B1 Light

0921-025


\$219.00 kit



Aegis® Orthodontic Adhesive with ACP




A light-cured orthodontic adhesive containing Amorphous Calcium Phosphate (ACP) that is used for bonding orthodontic brackets and bands to etched enamel. When orthodontic braces are removed, there may be signs of tooth demineralization in the area under or around the bracket. ACP may help prevent the demineralization that occurs due to micro-leakage or non-compliant patient oral hygiene. Aegis’ controlled placement keeps adhesive on the tooth while the practitioner secures accurate position of the bracket. Newly formulated material makes brackets stick perfectly!

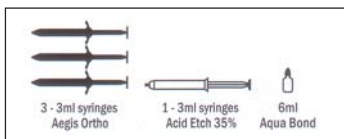

Standard Kit

0921-016

\$187.00 kit

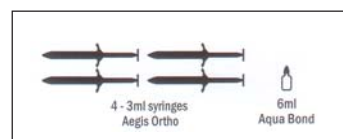

Refill Kit

0921-017


\$219.00 kit



Aegis® Pit & Fissure Sealant with ACP




The ACP filler is resilient and flexible, creating a stronger and longer lasting sealant. Aegis P&F Sealant has low shrinkage, good margin, and a slow rate of stress development so as not to chip off! Aegis P&F Sealant can be used with any light curing device and offers the added convenience of dispensing through syringe tips. The controlled flowability keeps the sealant on the tooth while completely filling the pits and fissures. When properly maintained and placed, Aegis forms a chemical and thermal barrier that can help keep patients free from carious lesions on the occlusal surfaces. Sealant cures to an opaque film on enamel. Standard kit contains enough material for approximately 30 patients.

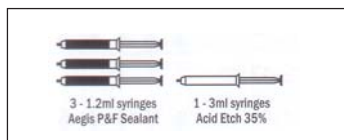

Standard Kit

0921-348

\$114.00 kit

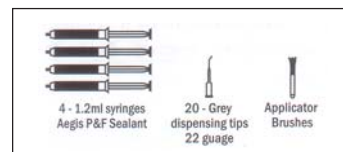

Refill Kit

0921-328

\$129.00 kit